



CSB1021HF LEC0131

FUNDAMENTALS OF GENOMIC DATA SCIENCE

0.0.0 Module 2: Assembly and annotation with Galaxy

0.1.0 About Fundamentals of Genomic Data Science

Fundamentals of Genomic Data Science is brought to you by the **Centre for the Analysis of Genome Evolution & Function (CAGEF)** bioinformatics training initiative. This course was developed based on feedback on the needs and interests of the Department of Cell & Systems Biology and the Department of Ecology and Evolutionary Biology.

The structure of this course is a “code-along”, hands-on style! A few hours prior to each lecture, materials will be made available for download at QUERCUS (<https://q.utoronto.ca/>). The teaching materials will consist of a weekly PDF that you can use to follow along with the instructor along with any datasets that you’ll need to complete the module. This learning approach will allow you to spend the time coding and not taking notes!

As we go along, there will be some in-class challenge questions for you to solve. Post lecture assessments will also be available for each module, building upon the concepts learned in class (see syllabus for grading scheme and percentages of the final mark).

0.1.1 Where is this course going?

We’ll take a blank slate approach here to learning genomic data science and assume you know nothing about programming or working directly with next generation sequencing data. From the beginning of this course to the end we want to guide you from potential scenarios like:

- You don’t know what to do with a set of raw sequencing files fresh from a facility like.

- You've been handed a legacy pipeline to analyse your data or maintain for the lab but you don't know what it runs or how.
- You plan on generating high-throughput data but there are no bioinformaticians around to help you out.

and get you to the point where you can:

- Recognize the basic tools in sequence analysis.
- Plan and write your own data analysis pipelines.
- Explain your data analysis methods to labmates, supervisors, and other colleagues.

0.1.2 How do we get there?

In the first half of this course, we'll focus on how to generate analysis pipelines using the Galaxy platform – a user-friendly graphical interface that provides access to common sequence analysis tools. After we are comfortable with these tools, we'll look at life through the lens of a command-line interface. It is here that we will learn the basics of file manipulation and how to program scripts that can carry out multiple tasks for us. From there we'll revisit tools from the first half and learn skills to make your data analysis life easier.

0.2.0 Goals of the module

1. Gain familiarity with the Galaxy interface.
2. Learn how to quality control sequencing reads with Galaxy.
3. Learn how to perform genome assembly and annotation with Galaxy.
4. Learn how to save and extract data and workflows with Galaxy.

0.3.0 Pre-class modules with Coursera

Each week we strongly encourage you to complete the assigned Coursera modules and/or readings **before** class. These are meant to provide you with sufficient background material on each week's module so that we can focus on the act of "doing" something with that data rather than spend a lot of time on the origins of it. You'll find a section outlining the next set of Coursera modules and readings at the end of each module.

0.3.1 Go to www.coursera.org and sign up for an account with your e-mail.

0.3.2 Search the following courses and enroll to audit each course (audit):

- Genomic Data Science with Galaxy, Johns Hopkins University.
- Command Line Tools for Genomic Data Science, Johns Hopkins University.

0.4.0 Setting up your working directory

We suggest that you create a new directory (folder) for this course directly off your root directory called **"FGDS"**. Working from your root directory is not necessary, but it will make some of the aspects of the course a little easier to manage. For MacOS users, we suggest you create this as a subfolder in your **user** directory.

- 0.4.1 Within this directory, create another directory called **"Module2"**. This is where we will store the data used in this week's module.
- 0.4.2 Create a subdirectory called **"downloads"** to store the initial files as we download them before decompressing and working with them in later steps.

1.0.0 Logging into Galaxy

This week we'll be working with the Galaxy platform and becoming familiar with how to navigate around it. There are a few ways you can use Galaxy between using an account over at <https://usegalaxy.org>, setting up a server instance using <https://genap.ca> or trying to set up your own Amazon Web services. Today we'll mainly cover the second option by working with an instance on GenAP.

Why are we using Galaxy? As we mentioned in the class slides, just as with the world of biological science, it's important that our experiments and results are reproducible. While we often encounter detailed methods in manuscripts, we might forget that the same applies to our data analysis.

While we may take this idea for granted, it's important to remember that software tools and packages can undergo revision! Running the same set of data between *different versions* of tools can result in separate outputs. Furthermore, if you don't list your parameter choices in detail, it can make your results harder to interpret and replicate.

While tracking all these details can be quite important, approaching software tools can also be daunting! We use Galaxy as a good way to ease into these software packages, learn to manipulate their parameters, and simplify the process of connecting their output to other tools. The graphical user interface greatly decreases the energy barrier to learning these common bioinformatic tools.


As we'll see later in class, you can also save these steps into workflows (pipelines) which you can use to repeat your analyses on different data sets or to tweak the parameters for a single dataset.

1.1.0 Galaxy Registration

After working with Galaxy in this course you may wish to explore it more using the publicly available Galaxy server architecture at <https://usegalaxy.org>. You'll find a number of tutorials, workflows, and datasets available here. After we're done with our work in this course you can use this as an alternate workspace to carry out some of your own pipelines.

1.1.1 Go to <https://usegalaxy.org>.

1.1.2 Select "Login or Register" located at the top of the page. Fill in all required fields and then submit.

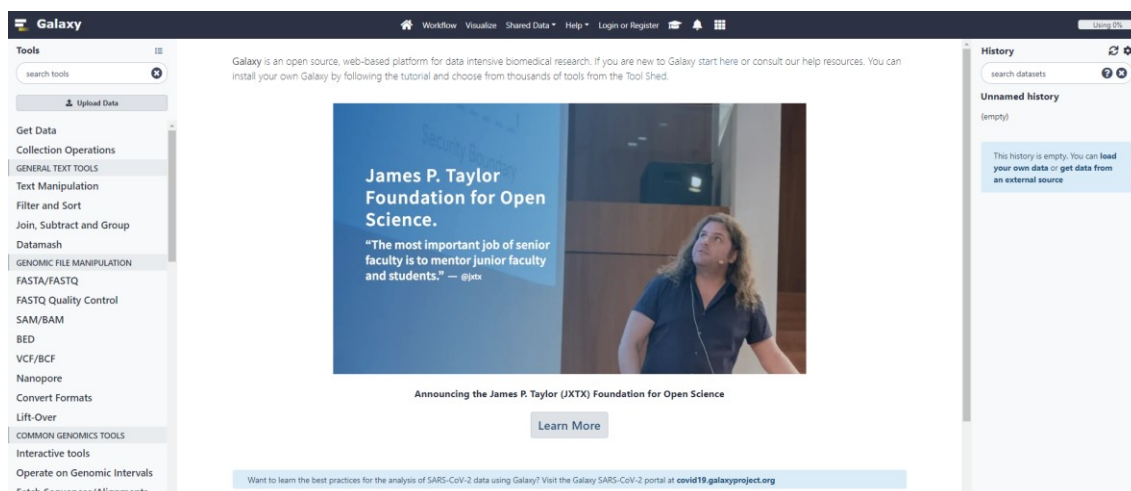


1.1.3 Verify your account by selecting the link in your email after it arrives.

The downside to using the online platform is that it can sometimes take several hours or even days for your analysis to get to the front of the queue


1.2.0 Running on our own galaxy instance through GenAP

GenAP is the Genetics & Genomics Analysis Platform and offers computing infrastructure to Canadian researchers. It is a cloud-based platform service that is built on Compute Canada resources. Therefore, having a Compute Canada Role Identifier (CCRI) can grant you easy access to these resources. For this course we have set up a galaxy instance at <https://galaxy-547-p38.p.genap.ca> and invited you to access it using your provided CCRI or a valid gmail address.

1.2.1 Login using your Compute Canada or Google credentials

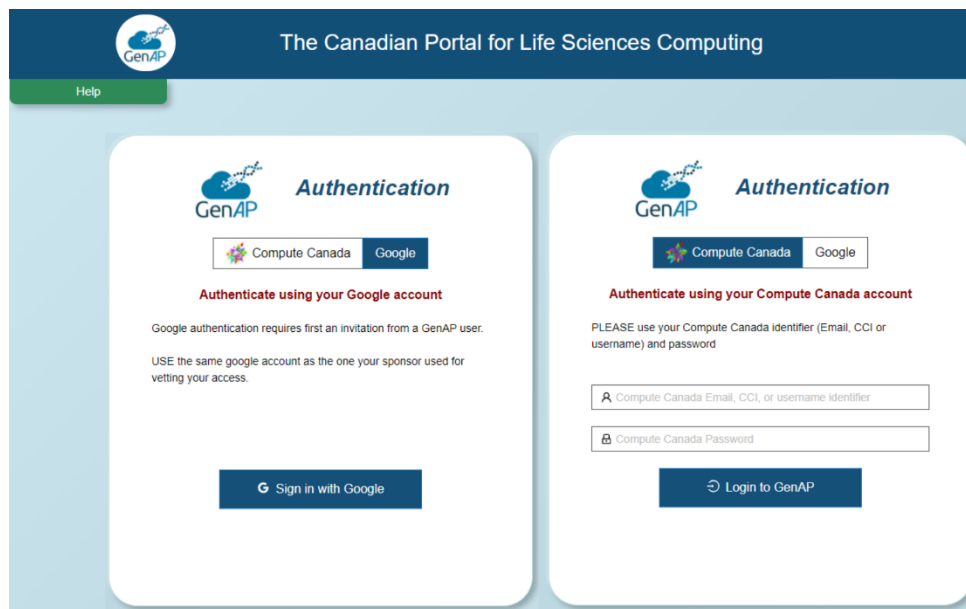

Once you are logged in, you'll be brought to the Galaxy homepage. It is from here that we'll be doing all of our activities.

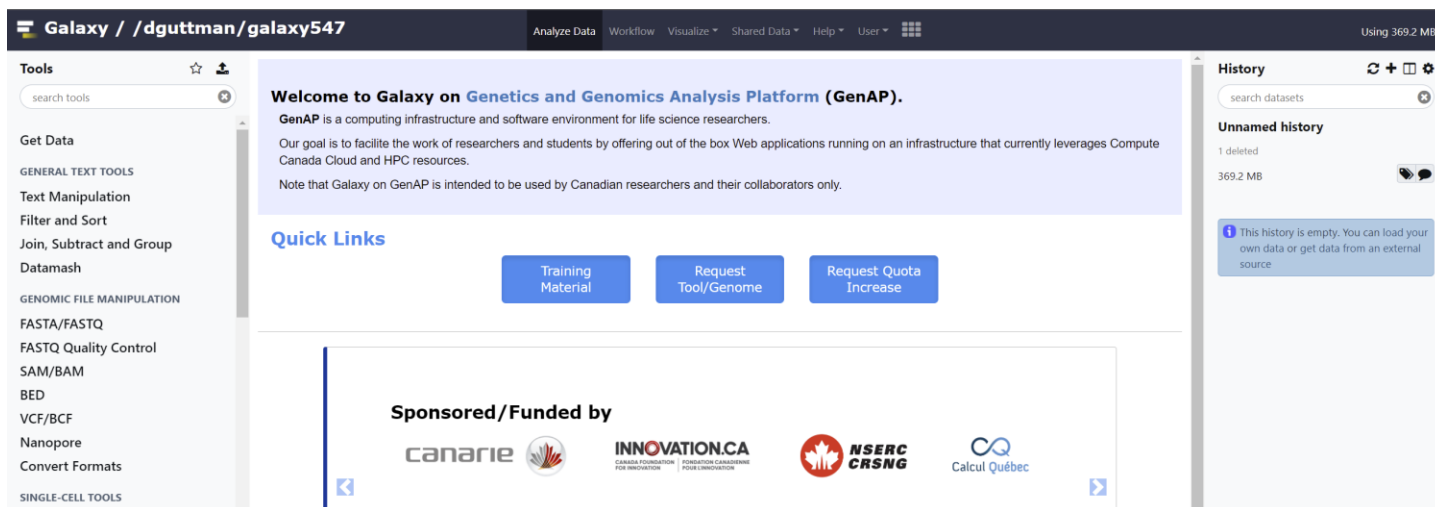

1.3.0 The Galaxy interface has four main quadrants

As you look at your galaxy home screen there are four quadrants to identify. Knowing the location and meaning of these will make navigation and following lecture instructions a little easier.

- 1.3.1 The Galaxy application screen is the main portion of your interface. It is centred on the screen and takes up 80% of the application interface. Here you can view results, datasets, and tool options as we navigate through our pipelines.
- 1.3.2 The Galaxy menu is located at the top of the screen. This menu can take you to a few different application screens where you can choose to Analyze Data, choose Workflows, and load Shared Datasets, as well as change user settings etc.
- 1.3.3 The Galaxy tools pane is located on the left side of the screen. This is most useful when working on the **Analyze Data** application screen. Scrolling down this tool menu you'll find popular tools. You can also use the Search bar to find tools by name.
- 1.3.4 The Galaxy History pane shows us the order of tools and applications that we have used. This can be used to see what you've done but also to put a pipeline together! We'll be referring to this often throughout our instructions.

Now that we've familiarized ourselves with the general interface of Galaxy, we can attempt to use some of the more popular tools and workflows. First, however, we need to import some data.



2.0.0 Importing datafiles in Galaxy

Before we can begin working with any data, we actually need to import it from somewhere. We can choose to upload data from local sources, from outside servers, or from within the Galaxy server itself. We'll walk through some helpful habits and instructions that will make it easier to import data for your data analysis

2.1.0 Creating a new history

- 2.1.1 In the History pane, choose the **+** symbol in the top right to create a new history.
- 2.1.2 Rename the history (directly below the history search bar) from **Unnamed history** to **AssemblyandAnnotation** and hit the Enter key.


2.2.0 Retrieving local data to analyse on Galaxy

For a typical analysis on Galaxy, you will upload the required data to the server using the **Get Data** tool from the tool pane. This option should be the top tool in your tool pane and clicking on it you'll see many options! This includes downloading and extracting FASTQ reads directly from the NCBI SRA but we won't do that today.

- 2.2.1 For now, we'll focus on **Get Data > Upload File from your computer**. Clicking on this option will bring up the following window:



- 2.2.2 The data you will need for this module is the raw fastq reads for *H. influenzae* str. KW20. Unfortunately, upload and download speeds can be painstakingly slow for large files. Therefore, we will instead do a mock upload with a much smaller file (**HinfKW20_rna.fna**) to provide an example of how this would work. Use the **Choose local file** button to upload this file data from last week.
- 2.2.3 Hit the **Start** button to begin the upload. Note that you can queue multiple files for upload!
- 2.2.4 While the file is uploading, you'll see it in the History pane first as a grey entry with a small clock icon, then as a yellowish entry with a spinning circle icon. When the operation is complete, the entry will turn a green colour.

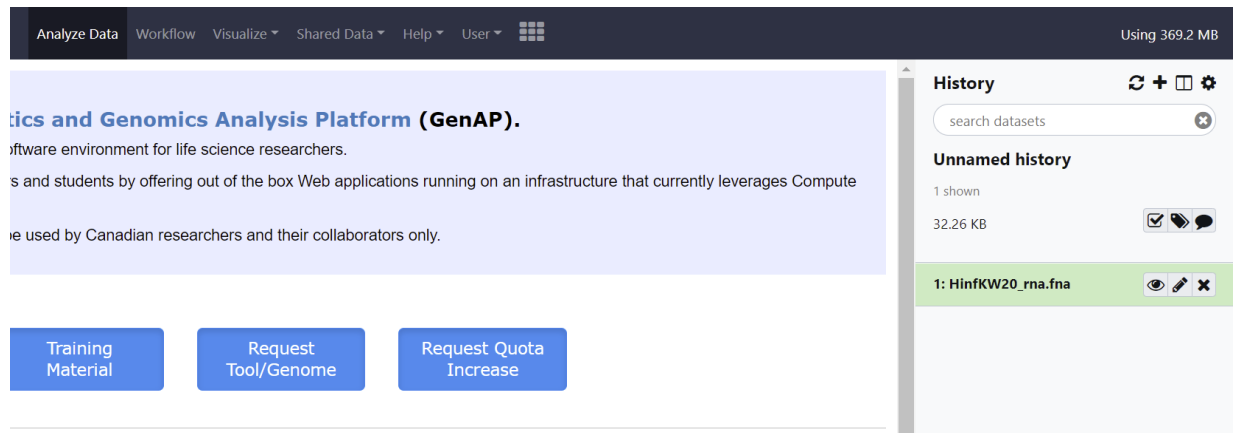


Notice that the entry has a number beside it? Each entry name in your *History pane* will be formatted as **#:history_entry** so that you can keep track of the order for your operations as well. As you create more entries, they will be listed in *descending* order.

2.3.0 Retrieving shared histories to analyse on Galaxy

While the above instructions are a useful example, we will still need the fastq file with our sequencing reads to do assembly and annotation. A subset file of only 2,000,000 reads has been uploaded to our server for you to access and load it to your own history using the *Shared Data* tab:

2.3.1. From the Galaxy *menu* choose **Shared Data > Histories**.

2.3.2 From the dataset list, choose **HinfKW20 fastq subset**. This will bring up the following information about the shared history:

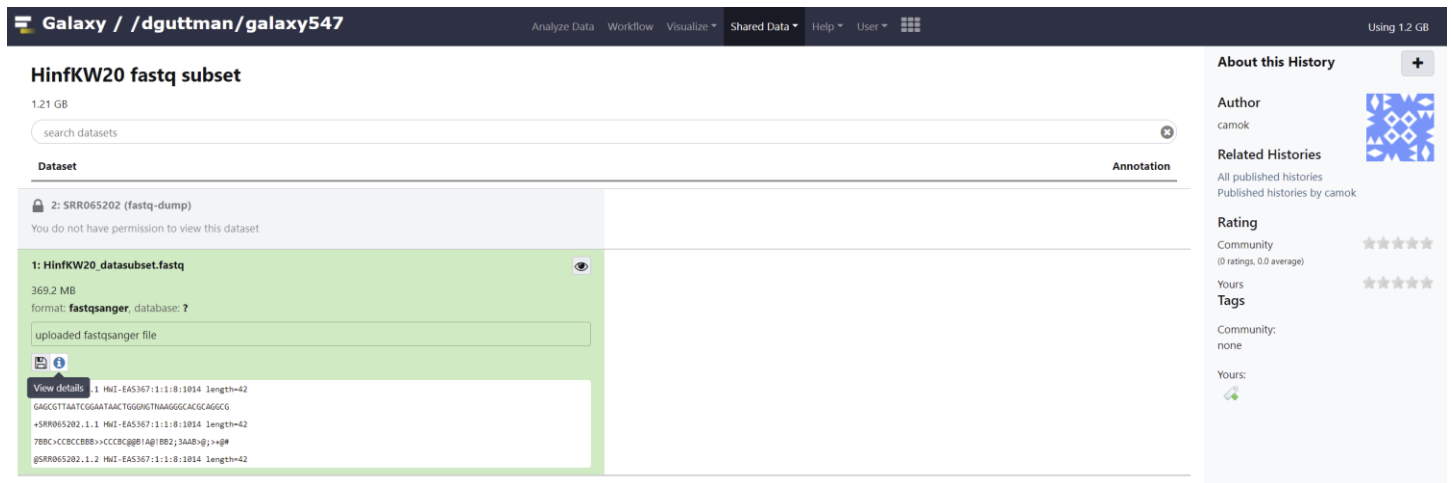

2.3.3 In the upper right corner you'll click on the large **+** to import the history to your Galaxy instance. This will bring the data into your version of Galaxy without creating multiple permanent copies across the class.


2.3.4 Name this imported history as **HinfKW20 fastq subset**. This action will return you to the main Galaxy homepage.

The process for adding shared datasets to your history will be similar. However, the current setup for the Galaxy instance facilitates shared data through *Shared Data > Histories* instead.

2.4.0 Add an imported history to your own history


Now that we've imported a history into our instance of Galaxy we want to add that to our current **AssemblyandAnnotation** history.

- 2.4.1 Click on the **View all histories** icon (double rectangle) at the top of the *History pane*. This will bring you to the history page which provides a fast way to switch between multiple histories. Here we can see the newly imported history as well as our **AssemblyandAnnotation** history.

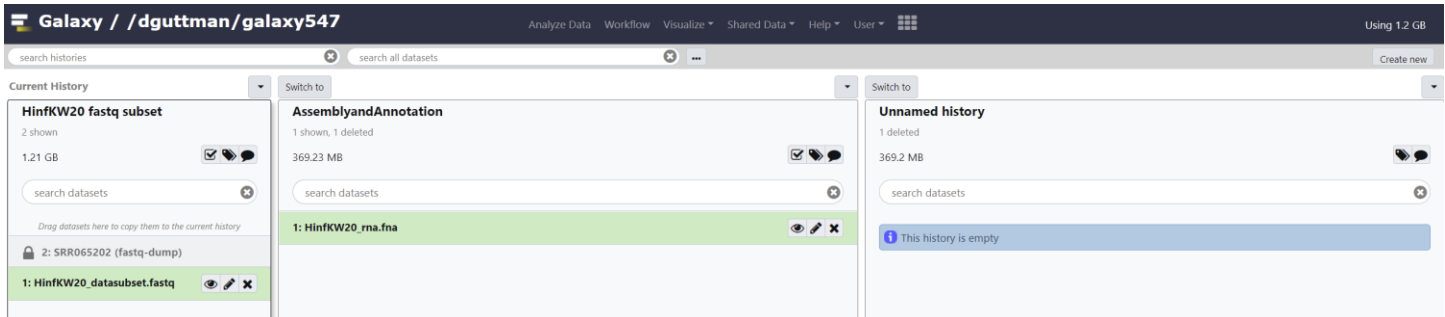

- 2.4.2 Above the **AssemblyandAnnotation** history, click on the **Switch to** button.

- 2.4.3 In the **HinfKW20 fastq subset**, select the **HinfKW20_datasubset.fastq** and add it to the **AssemblyandAnnotation** history by dragging and dropping it into the supplied box.

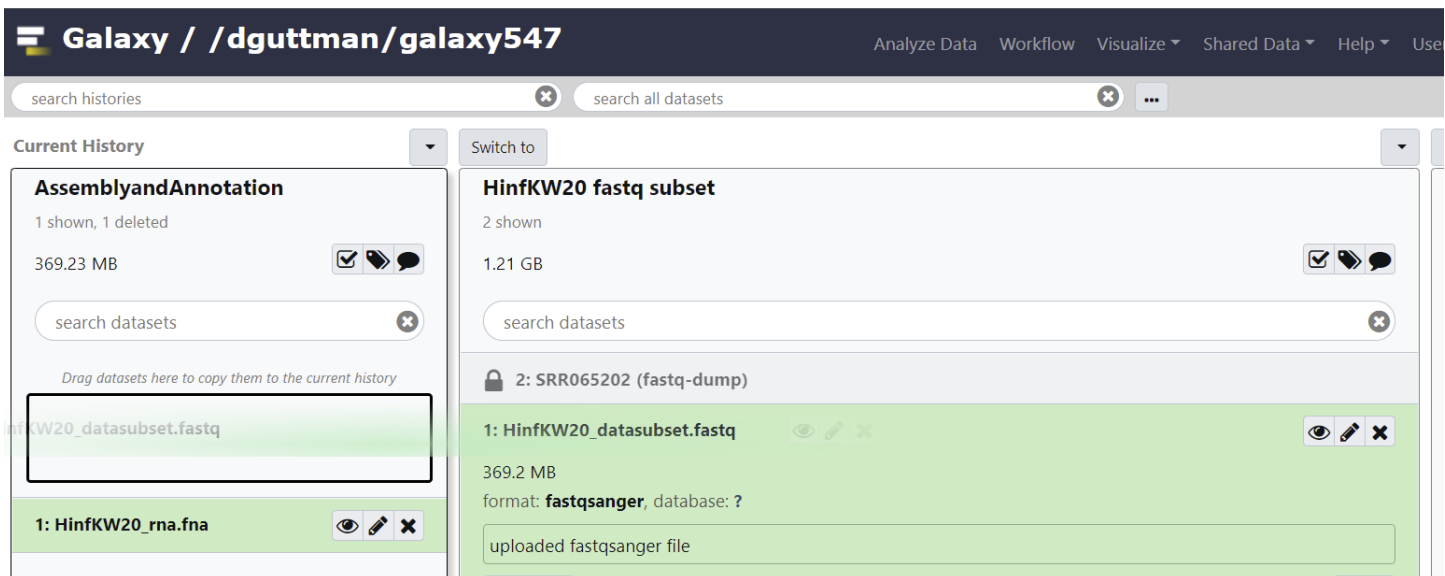

- 2.4.4 Return to the Galaxy homepage by clicking on the top left corner.

- 2.4.5 You will now see the **HinfKW20_datasubset.fastq** file in your "History" panel. Verify that the datatype is **fastqsanger** using the Edit Attributes (pencil) key and clicking on the Datatypes tab.

- 2.4.6 Change datatype to **fastqsanger** if necessary.

3.0.0 Quality filtering raw sequencing reads

Now that our data subset is available in our history, the first step is to perform a quality inspection with **FastQC** to identify potential issues with our data. After quality inspection, we will perform quality control with **Trimmomatic** (if necessary) to do to get rid of the identified issues.

3.1.0 Run *FastQC* on raw sequencing reads

The FastQC tool is a popular and helpful way to check the overall read quality of your dataset. While some sequencing platforms like Illumina will perform a set of quality filtering as they run, these are really just eliminating reads with a low passing filter. This means, as a whole the reads show quality issues that make their overall quality questionable. The FastQC tool takes this a step further to look closely at the Q-scores in a number of modules to create a human-readable html report.

3.1.1 From the Galaxy [tools pane](#) select **FASTQ Quality Control > FastQC Read Quality reports**

3.1.2 Scroll down to review the purpose and input/output information for *FastQC*.

3.1.3 Execute *FastQC* with the following settings:

- Short read data from your current history: **HinfKW20_datasubset.fastq**
- Other options: **default**

The run will generate two new history entries:

- **3: FastQC on data 2: Webpage**
- **4: FastQC on data 2: RawData**

Note that in some versions, you may be given the option to title your output. Our GenAP instance uses the history value (**data 2**) as the naming scheme to produce output.

The screenshot displays the Galaxy web interface. The top navigation bar includes 'Galaxy / /dguttman/galaxy547', 'Analyze Data', 'Workflow', 'Visualize', 'Shared Data', 'Help', and 'User'. The left sidebar shows a 'Tools' pane with a search bar and categories like 'Get Data', 'GENERAL TEXT TOOLS', 'Text Manipulation', 'Filter and Sort', 'Join, Subtract and Group', 'Datamash', 'GENOMIC FILE MANIPULATION', 'FASTA/FASTQ', and 'FASTQ Quality Control'. The main panel shows a green message: 'Executed FastQC and successfully added 1 job to the queue.' It lists the input as '2: HinfKW20_datasubset.fastq' and the outputs as '4: FastQC on data 2: RawData' and '3: FastQC on data 2: Webpage'. A note at the bottom states: 'You can check the status of queued jobs and view the resulting data by refreshing the History panel. When the job has been run the status will change from 'running' to 'finished' if completed successfully or 'error' if problems were encountered.' The right sidebar shows the 'History' pane with a search bar and a list of jobs. The top job is '4: FastQC on data 2: Raw Data' (369.23 MB). Below it is '3: FastQC on data 2: Webpage'. The bottom job is '2: HinfKW20_datasubset.fast' (369.23 MB). The bottom job is highlighted in green.

3.2.0 Reviewing your *FastQC* output

Once the background of the output files is green, the analysis is complete. We can review the results of the FastQC analysis on the [History pane](#).


3.2.1 Review the HTML output from **4: FastQC on data 2: Webpage** by selecting the [View data](#) icon (eye) on the history entry.

3.2.2 You will see the following result sections that you can review:


- Basic Statistics (Pass)
- Per base sequence quality (Pass)
- Per tile sequence quality (Pass)
- Per sequence quality scores (Pass)
- Per base sequence content (Warning)
- Per base GC content (Pass)
- Per sequence GC content (Pass)
- Per base N Content (Pass)
- Sequence Length Distribution (Pass)
- Sequence Duplication Levels (Warning)
- Overrepresented Sequences (Pass)
- Adapter Content (Pass)

3.2.3 From our analysis we can see that overall, the read quality of along the set of raw reads remains quite high throughout the length of the reads. Even towards the tail end of the reads, the median Q score remains above 28. The average read quality per sequence is ~32.5 and the distribution from 2M reads shows bulk of reads to have a mean Q score above 28.

✔ Per base sequence quality

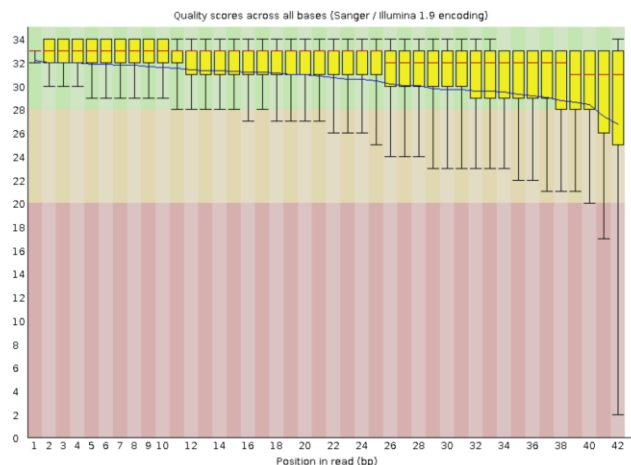


✔ Per sequence quality scores

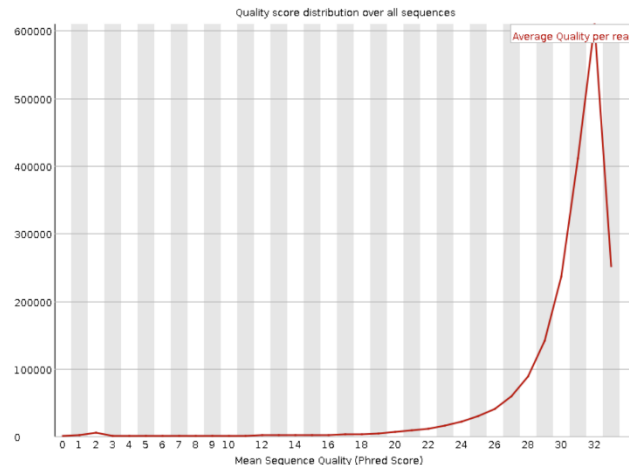


3.2.4 While this analysis reveals that there are no major problems with our dataset, although sequence content is slightly biased at the start of every read, and there is a slight excess of duplicated sequences:

! Per base sequence content

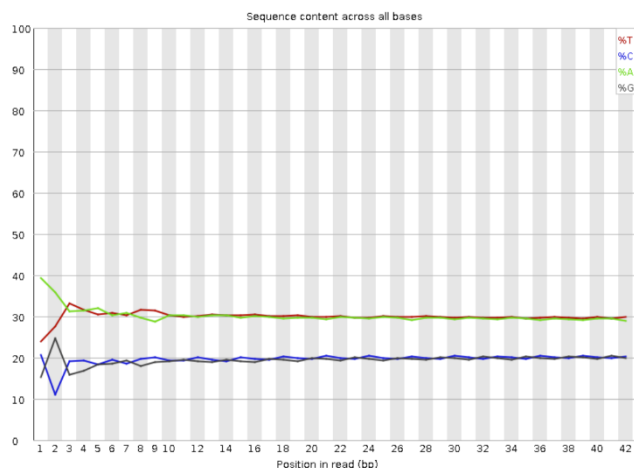

! Sequence Duplication Levels

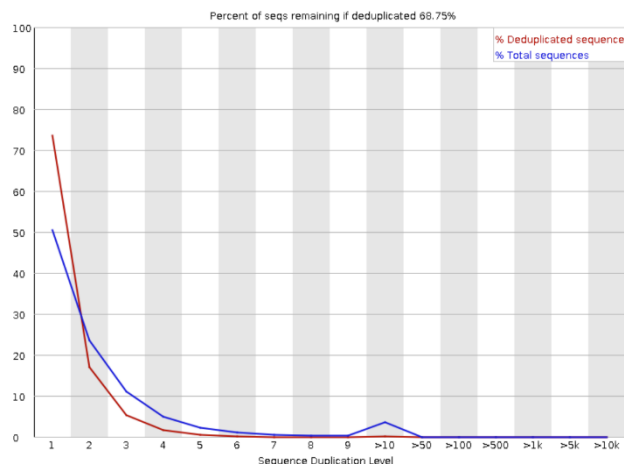

3.3.0 Use *Trimmomatic* to remove poor-quality reads

Trimmomatic performs a variety of useful tasks to help remove unwanted sequences from your reads. Developed for Illumina sequencing reads, this tool can remove adapter sequences generated during library preparation, while also removing sequences from the ends of your reads based on criteria like quality score.

Using this tool on our dataset may help remove some of the warnings received during our *FastQC* analysis. However, care should always be taken when trimming or removing reads, because strict filters can result in a lot of lost data and create new issues with sequencing reads.

3.3.1 From the Galaxy *tools pane* select **FASTQ Quality Control > *Trimmomatic* flexible read trimming tool for Illumina NGS data**.

3.3.2 Review the parameters that can be used for *Trimmomatic* and the required inputs. Some commonly employed parameters include:

- **ILLUMINACLIP:** Trim adaptor and other Illumina sequences off the ends of reads.
- **SLIDINGWINDOW:** Trim reads based on average quality in a window of defined size.
- **MINLEN:** Drop all reads below a defined length.
- **LEADING/TRAILING:** Cut bases at the start/end of reads if below defined quality.
- **CROP/HEADCROP:** Cut a defined number of bases at start/end of reads (regardless of quality).
- **AVGQUAL:** Drop all reads below a defined average quality.

3.3.3 Execute Trimmomatic with the following settings:

- Single-end or paired-end reads: **Single-end**
- Input FASTQ file: **2: HinfKW20_datasubset.fastq**
- Perform initial ILLUMINACLIP step? **Yes** (generates the following default values)
 - Select standard adapter sequences or provide custom? **Standard**
 - Adapter sequences to use: **TruSeq2 (paired-ended for Illumina GAII)**
 - Maximum mismatch count: **2**
 - Minimum matches between two 'adapter ligated' reads: **30**

- Minimum matches between adapter sequence and a read: **10**
 - Minimum adapter length that needs to be detected: **8**
- Trimmomatic Operation:
 - Sliding window trimming (SLIDINGWINDOW)** with default values
 - Number of bases to average across: **4**
 - Average quality required: **20**
- Insert Trimmomatic Operation:
 - Cut the specified number of bases from the start of the read (HEADCROP)**
 - Number of bases to remove from the start of the read: **10**

- 3.3.4 Note that with **Trimmomatic**, the order of operations matters. For example, if you performed the SLIDINGWINDOW trimming before using CROP, you would be cropping reads that have already been trimmed. This would not usually be what you want, because the SLIDINGWINDOW will already have removed many of the low-quality bases at the end of reads.
- 3.3.5 To view some of the run details, in the History pane, click on the entry for **Trimmomatic on HinfKW20_datasubset.fastq**. This will expand the history entry.
- 3.3.6 Click on the View details (circle with i) icon which will bring up the running parameters of the Trimmomatic run

```

5: Trimmomatic on HinfKW20
_datasubset.fastq
314.3 MB
format: fastqsanger, database: ?

TrimmomaticSE: Started with arguments:
-threads 1 -phred33 fastq_in.fastqsanger
fastq_out.fastqsanger
ILLUMINACLIP:/cvmfs/soft.galaxy/v2.1/tool-
dependency/_conda/envs/_trimmomatic@0.
0.36-6/adapters/TruSeq2-SE.fa:2:30:10
SLIDINGWIND

```

View details

```

.1 HWI-EAS367:1:1:8:1014 length=42
CGGAATAACTGGGNGT
+SRR065202.1.1 HWI-EAS367:1:1:8:1014 length=42
BBB>>CCCB@B!A@
+SRR065202.1.2 HWI-EAS367:1:1:8:1014 length=42

```

3.3.7 Under the Job Information section, find the Tool Standard Output line and select **stdout**.

Galaxy / /dguttman/galaxy547

Analyze Data Workflow Visualize Shared Data Help User

Tools

search tools

Get Data

GENERAL TEXT TOOLS

Text Manipulation

Filter and Sort

Join, Subtract and Group

Datamash

GENOMIC FILE MANIPULATION

FASTA/FASTQ

FASTQ Quality Control

MultiQC aggregate results from bioinformatics analyses into a single report

Trimmomatic

Dataset Information

Number: 6

Name: Trimmomatic on HinfKW20_datasubset.fastq

Created: Mon Nov 8 20:55:31 2021 (UTC)

Filesize: 314.3 MB

Dbkey: ?

Format: fastqsanger

Job Information

Galaxy Tool ID: toolshed.g2.bx.psu.edu/repos/pjbriggs/trimmomatic/trimmomatic/0.36.5

Galaxy Tool Version: 0.36.5

Tool Version:

Tool Standard Output: **stdout**

Tool Standard Error: stderr

Tool Exit Code: 0

History Content API ID: 2a56795cad3c7db3

Job API ID: a799d38679e985db

History API ID: ebf8f50c6abde6d

3.3.8 Note that of 2M reads, 96.94% survived while 3.06% were dropped. The **stdout** (standard out) information also includes the Illumina adapter sequences used in the clipping process.

3.4.0 Rerun your trimmed data in FastQC


Now that we've trimmed our data of adapters sequences and some poor-quality reads, we can return to examine how this changes our FastQC output.

3.4.1 From the Galaxy tools pane select **FASTQ Quality Control > FastQC Read Quality reports**

3.4.2 Execute FastQC with the following settings:


- Short read data from your current history: **Trimmomatic on HinfKW20_datasubset.fastq**
- Other options: **default**

3.4.3 Review the HTML output from **6: FastQC on data 5: Webpage** by selecting the View data icon (eye) on the history entry.

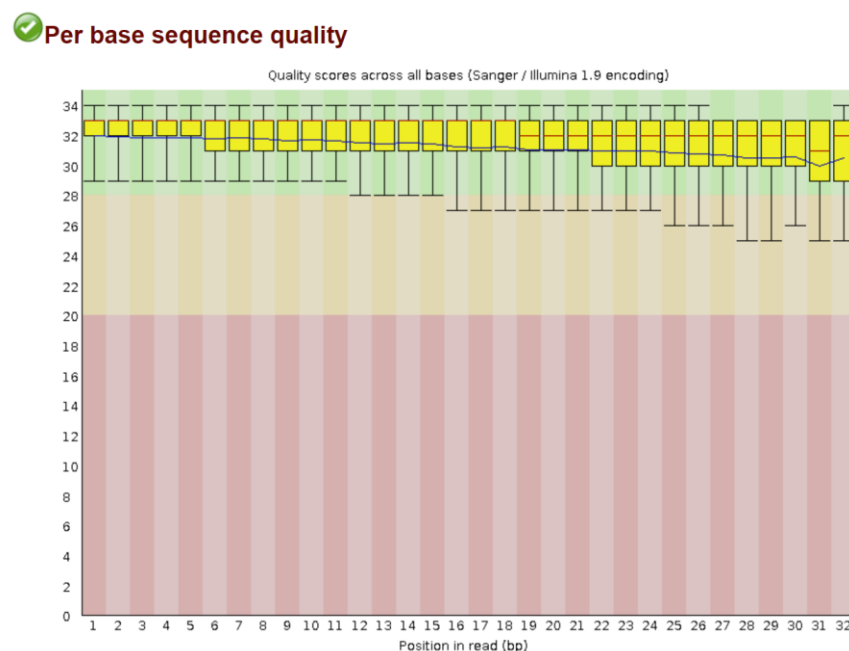


- We have solved our initial issues with **Per base sequence content** through the **HEADCROP** command which removed the first 10 bases of each read thus reducing the initial biases at the start of our reads.

✔ **Per base sequence content**

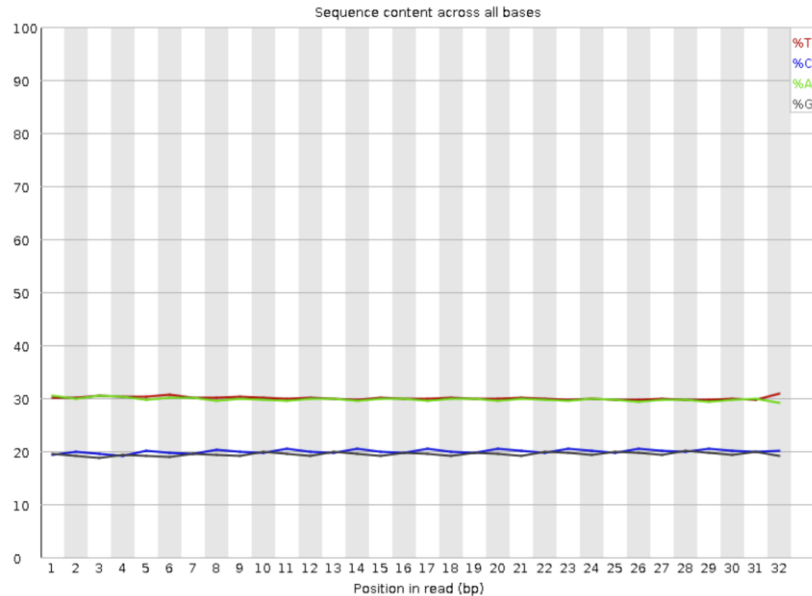

- As a consequence, however, our reads are now only 32 bases and we get a warning for the **Per sequence GC content** as the distribution is shifting slightly to a smaller standard deviation.

⚠ **Per sequence GC content**

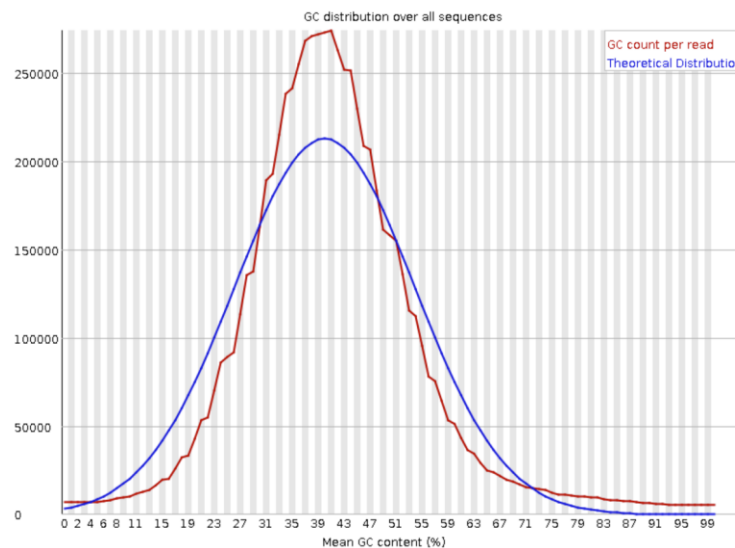

- We can also see from our **Sequence Length Distribution** that our reads are no longer strictly 42 basepairs but range from 14 to 32 basepairs in length.

4.0.0 Genome Assembly with Velvet

To recap, we have imported raw sequence reads from the *Haemophilus influenza* genome and have trimmed and filtered these reads using a combination of *Trimmomatic* and *FastQC* to analyse the read quality and distribution. We are now ready to assemble these reads into our first draft genome assembly with the Velvet tool.

Given the nature of our reads now ranging in size from 14 to 32 basepairs we can assume that we will end up with several smaller contigs (assembled fragments). Other contributors of small contigs can be low coverage or repetitive regions. Regardless, we'll have an assembly that we can begin to annotate in later steps.

The *Velvet* assembly tool runs in two steps by

- 1) Use *velveth* to build a **h**ash table of k-mers in each read ie. it tracks each read by its different subsequences.
- 2) Use *velvetg* to perform a **g**enome assembly based on the hash table from *velveth*

4.1.0 Use *velveth* to produce a hash table of sequences

- 4.1.1 From the Galaxy *tools pane* select **Assembly > *velveth* Prepare a dataset for the Velvet velvetg Assembler**.
- 4.1.2 Scroll down to review the purpose and input/output information about *velveth*.
- 4.1.3 Execute *velveth* to create your hash table using the following settings:
 - Hash Length: **21**
 - Use strand specific transcriptome sequencing: **No**
 - Input Files:
 - Choose the input type: **Single ended**
 - Read type: **short reads**
 - Dataset: **Trimmomatic on HinfKW20_datasubset.fastq**

The screenshot shows the Galaxy web interface for the tool 'velveth Prepare a dataset for the Velvet velvetg Assembler (Galaxy Version 1.2.10.3)'. The interface is divided into a top navigation bar, a left-hand 'Tools' sidebar, and a main configuration area. The 'Tools' sidebar lists various tools, including 'velveth Prepare a dataset for the Velvet velvetg Assembler'. The main configuration area has the following settings:

- Hash Length:** 21 (with a slider and a note: 'k-mer length in base pairs of the words being hashed.')
- Use strand specific transcriptome sequencing:** No (with a note: 'If you are using a strand specific transcriptome sequencing protocol, you may wish to use this option for better results.')
- Input Files:** 1: Input Files
- Choose the input type:** Single ended
- read type:** short reads
- Dataset:** 6: Trimmomatic on HinfKW20_datasubset.fastq

4.2.0 Use *velvetg* to produce a genome assembly

Now that we have our hash table from *velveth* we can send this on the *velvetg* which will use the information to construct a de Bruijn graph by adding the k-mers one-by-one to the graph and then chaining and merging paths and nodes to produce the final contig sequences.

- 4.2.1 From the Galaxy [tools pane](#) select **Assembly > velvetg Velvet sequence assembler for very short reads**.
- 4.2.2 Scroll down to review the purpose and input/output information about [velvetg](#).
- 4.2.3 Execute [velvetg](#) to produce your assembly using the following settings:
 - Velvet Dataset: **velveth on data 5**
 - Coverage cutoff: **Automatically determined**
 - Tracking of short read positions in assembly: **No**
 - Set minimum contig length: **Yes**
 - Minimum contig length: **200**
 - Expected Coverage of Unique Regions: **None**
 - Additional outputs: **check this box** (if you don't you'll get a Java error)
 - Generate a AMOS.afg file
 - Generate a UnusedReads fasta file
 - Generate a velvet LastGraph file
 - Using Paired Reads: **No**

Galaxy / /dguttman/galaxy547

Analyze Data Workflow Visualize Shared Data Help User

Tools ☆

search tools

scRNA-seq: Scanpy (IUC)
scRNA-seq: STARsolo
scRNA-seq: UCSC Cell Browser

COMMON GENOMICS TOOLS
Operate on Genomic Intervals
Fetch Sequences/Alignments

GENOMICS ANALYSIS
Assembly
[velvetg](#) Velvet sequence assembler for very short reads
[velveth](#) Prepare a dataset for the Velvet velvetg Assembler
[Create assemblies with Unicycler](#)
[Quast](#) Genome assembly Quality

Annotation
Mapping
Variant Calling
ChIP-seq
RNA-seq
Multiple Alignments
Phenotype Association
Regional Variation
STR-FM: Microsatellite Analysis

<https://galaxy-547-p38.p.genap.ca>

velvetg Velvet sequence assembler for very short reads (Galaxy Version 1.2.10.2) ☆ Favorite Options

Velvet Dataset
8: velveth on data 5

Prepared by velvetg.

Coverage cutoff
None

Tracking of short read positions in assembly
Yes No

Generates Graph2 dataset

Set minimum contig length
Yes
minimum contig length exported to contigs.fa file (default: hash length * 2).
minimum contig length
200
minimum contig length exported to contigs.fa file (default: hash length * 2)

Expected Coverage of Unique Regions
None

Additional outputs
☒ Select/Unselect all
☒ Generate a AMOS.afg file ☒ Generate a UnusedReads fasta file ☒ Generate velvet LastGraph file

Using Paired Reads
No

4.3.0 Reviewing your velvetg results

Looking at our [History pane](#), we see that [velvetg](#) has generated five sets of output but we are just interested in two: **velvetg on data 8: Contigs** and **velvetg on data 8: Stats**.

- 4.3.1 Review the **stdout** file from [velvetg on data 8: Contigs > View Details > stdout](#).

You will be able to find an N50 and total length for your assembly at the bottom of file, but these values are for your de Bruijn graph so don't account for contig filtering. N50 is a quality metric that is defined as

the shortest sequence length at 50% of the genome. Higher N50 values mean that you have a more contiguous assembly.

4.3.2 Review the **Stats** output from velvetg on data 8: Stats > View Data.

| ID | lgth | out | in | long_cov | short1_cov | short1_Ocov | short2_cov | short2_Ocov | short3_cov | short3_Ocov | short4_cov | short4_Ocov | long_nb | short1_nb | short2_nb | sl |
|----|------|-----|----|----------|------------|-------------|------------|-------------|------------|-------------|------------|-------------|---------|-----------|-----------|----|
| 1 | 2111 | 0 | 0 | 0.000000 | 14.331596 | 14.331596 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 2 | 1923 | 0 | 0 | 0.000000 | 12.835153 | 12.835153 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 3 | 1456 | 0 | 0 | 0.000000 | 11.760989 | 11.760989 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 4 | 2052 | 0 | 1 | 0.000000 | 17.330897 | 17.330897 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 5 | 3773 | 1 | 0 | 0.000000 | 13.667638 | 13.667638 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 6 | 352 | 0 | 0 | 0.000000 | 6.519886 | 6.519886 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 7 | 1562 | 0 | 0 | 0.000000 | 19.039052 | 19.039052 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 8 | 1126 | 0 | 0 | 0.000000 | 9.671403 | 9.671403 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 9 | 2100 | 2 | 1 | 0.000000 | 13.714762 | 13.714762 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 10 | 2240 | 2 | 0 | 0.000000 | 19.016071 | 19.016071 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 11 | 504 | 0 | 0 | 0.000000 | 7.380952 | 7.380952 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 12 | 4955 | 0 | 0 | 0.000000 | 13.040565 | 13.040565 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 13 | 1151 | 0 | 1 | 0.000000 | 9.097307 | 9.097307 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 14 | 901 | 0 | 0 | 0.000000 | 11.073252 | 11.073252 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 15 | 1000 | 1 | 0 | 0.000000 | 16.652000 | 16.652000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 16 | 1105 | 0 | 0 | 0.000000 | 9.086878 | 9.086878 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 17 | 1298 | 0 | 0 | 0.000000 | 16.143297 | 16.143297 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 18 | 1627 | 1 | 0 | 0.000000 | 17.536570 | 17.536570 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 19 | 148 | 0 | 2 | 0.000000 | 18.898649 | 17.371622 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 20 | 4338 | 0 | 1 | 0.000000 | 14.173352 | 14.173352 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 21 | 3932 | 1 | 1 | 0.000000 | 25.207528 | 25.207528 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 22 | 3183 | 1 | 0 | 0.000000 | 13.555765 | 13.555765 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 23 | 705 | 0 | 0 | 0.000000 | 9.235461 | 9.235461 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 24 | 1946 | 0 | 0 | 0.000000 | 12.951182 | 12.951182 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 25 | 358 | 0 | 0 | 0.000000 | 6.536313 | 6.536313 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |
| 26 | 813 | 0 | 1 | 0.000000 | 26.521525 | 26.521525 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0.000000 | 0 | 0 | 0 | 0 |

Now that we have our genome assembled, and a list of contigs, what should we do with it? Annotate it of course!

What is a de Bruijn graph? Yes, you have generated a genome using Velvetg but how does it *actually* work? All the **k-mers** are represented in a graph with connected neighbours being **k-mers** overlapping by **(k-1)** bases. Each node stores information about how often it's k-mer is seen within the data. For more information on using the Velvet *de novo* assembler check out [Zerbino 2010 on Pubmed](#)

5.0.0 Genome annotation

As is the case with most forms of genome analysis, there are a large number of options available for genome annotation - the process of identifying coding regions in a genome and assigning a function to them.

5.1.0 Genome annotation with Prokka

Prokka is designed for smaller genomes like those of bacteria, archaea, and viruses. Its major advantages are speed and simplicity. **Prokka** can fully annotate a typical 4Mbp genome in less than 10 minutes on most modern personal computers. Furthermore, it integrates coding sequence identification and functional assignment into a single command, which makes it very easy to use.

For eukaryotic genomes, you must consider introns and exons for each gene. It is also important that you mask repetitive sequences in these genomes. Some common programs for the annotation of eukaryotic genomes include **Augustus**, **Maker**, and **Pasa**. Information about these programs is available online.

5.1.1 From the Galaxy [*tools pane*](#) select **Annotation > Prokka Prokaryotic genome annotation**

5.1.2 Scroll down to review the purpose and input/output information about Prokka.

There are lots of options here that can be useful depending on the goals of the annotation. These include, for example, specifying genus or species-specific databases to prioritize when assigning functions to coding regions. You can also specify which output files you want Prokka to generate. Today, we will simply be running Prokka using default settings and will be assigning functions with the default BLAST database.

5.1.3 Execute Prokka to produce your annotation using the following settings:

- Contigs to annotate: **velvetg on data 8: Contigs**
- Locus tag prefix: **HinfKW20Gene**
- Genus name: **Haemophilus**
- Species name: **influenzae**
- Strain name: **KW20**

5.1.4 Scrolling past the remaining parameters (leave as **default**) we find the [*Additional outputs*](#) section where we can choose which files to create from the annotation. Leave these as the default (**Select all**).

The screenshot shows the Galaxy web interface. At the top, the header indicates the user is logged in as 'dguttman/galaxy547'. The main content area displays a green message box stating 'Executed Prokka and successfully added 1 job to the queue.' Below this, it lists the input: '13: velvetg on data 8: Contigs' and the outputs: '25: Prokka on data 13: log' and '24: Prokka on data 13: txt'. On the left, the 'Tools' pane is visible, showing a search bar and a list of tools under 'COMMON GENOMICS TOOLS'. On the right, the 'History' pane shows a list of datasets, including '25: Prokka on data 1 3: log' and '24: Prokka on data 1 3: txt'.

5.2.0 Review your Prokka output

You'll see that our **Prokka** execution has generated 12 files, many of which should appear familiar to you from module 1: gff, gbk, fna, faa, ffn, sqn, fsa, tbl. In addition to this it has also produced a coding sequence (tsv); a discrepancy report (err); a brief summary (txt); and a more verbose run report (log).

Many of these files are what we would find accompanying our assembly report from NCBI.

6.0.0 Downloading data and saving your pipeline


Downloading output files and applying appropriate naming conventions from Galaxy analyses is not as easy as we might like. This process can be sped up somewhat by gaining familiarity with the command line, but in this section of the course we will stick with **Finder** and **File Explorer**.


6.1.0 Create a dataset list of your data from the *History* pane


- 6.1.1 Recall that you have already set up a subfolder for this module under your main FGDS folder.
- 6.1.2 In the Galaxy History pane select the Operations on multiple datasets (checkbox) button. You'll see selectable checkboxes are now available for each history entry.
- 6.1.3 Select the following history entries to download:



- 6.1.4 From the *History pane* select the For all selected... dropdown menu and select **Build Dataset List**

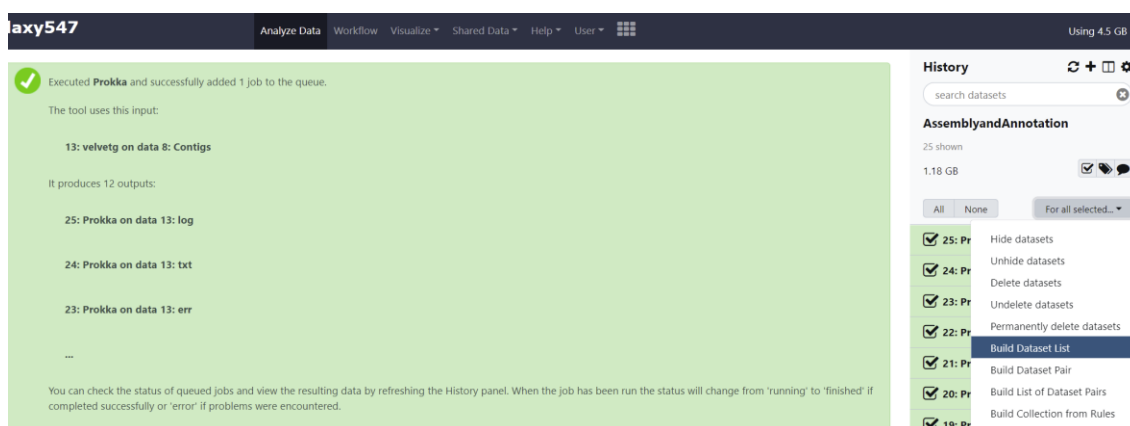

- 6.1.5 Name the list  You will now see a new entry in the *History pane* with the name of your newly created list.

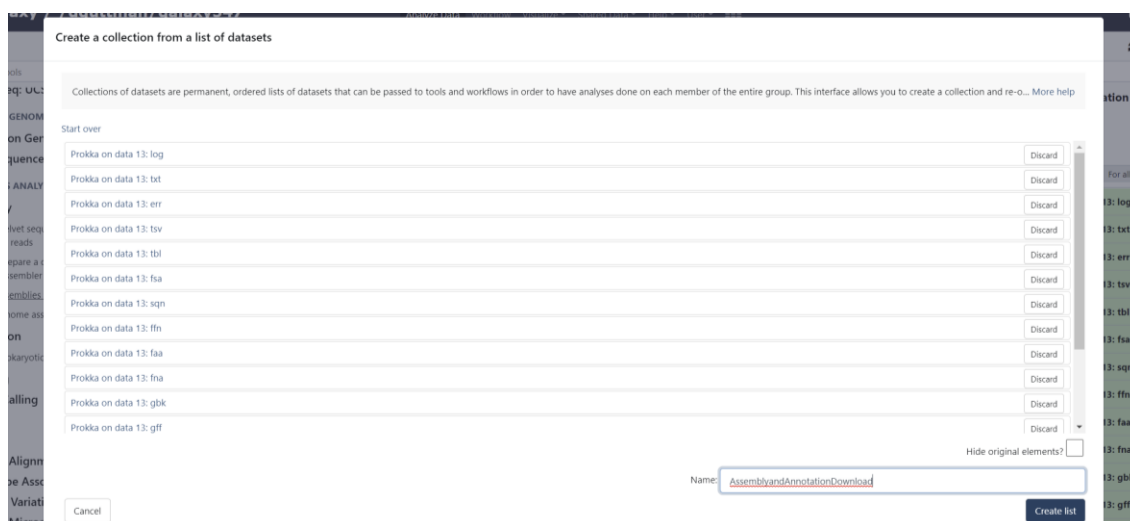

6.2.0 Download your dataset list

Now that we've created a list of 15 items we can download it to your own desktop.

- 6.2.1 Go to the *History pane* and deselect the *Operations on multiple datasets* checkbox icon (if needed) and then select the **AssemblyandAnnotationDownload** entry. This will show the 15 items of your list.
- 6.2.2 From the *History pane* click on the *Download Collection* button and save the file to your **~/FGDS/Module2/downloads** folder. Notice that a single **.tgz** archive file will download to your directory.
- 6.2.3 Decompress the **AssemblyandAnnotationDownload.tgz** into **AssemblyandAnnotationDownload.tar** in your favourite program.
- 6.2.4 Unarchive **AssemblyandAnnotationDownload.tar** and place its contents into **~/FGDS/Module2/** and rename the files (and folders) as follows:

| Download name | New name |
|----------------------------------|------------------------------------|
| FastQC on data 2_ Webpage.html | FastQC_RawOutput.html |
| FastQC on data 2_ Webpage | FastQC_RawOutput |
| FastQC on data 4_ Webpage.html | FastQC_Trimmed.html |
| FastQC on data 4_ Webpage | FastQC_Trimmed |
| Prokka on data x_ err.txt | HinfKW20_prokka.err |
| Prokka on data x_ faa.fasta | HinfKW20_prokka.faa |
| Prokka on data x_ ffn.fasta | HinfKW20_prokka.ffn |
| Prokka on data x_ fna.fasta | HinfKW20_prokka.fna |
| Prokka on data x_ fsa.fasta | HinfKW20_prokka.fsa |
| Prokka on data x_ gbk.txt | HinfKW20_prokka.gbk |
| Prokka on data x_ gff.gff | HinfKW20_prokka.gff |
| Prokka on data x_ log.txt | HinfKW20_prokka.log |
| Prokka on data x_ sqn.asn1 | HinfKW20_prokka.sqn.asn1 |
| Prokka on data x_ tbl.txt | HinfKW20_prokka.tbl |
| Prokka on data x_ tsv.tabular | HinfKW20_prokka_cds.tsv |
| Prokka on data x_ txt.txt | HinfKW20_prokka.txt |
| velvetg on data x_ Contigs.fasta | HinfKW20_denovo_assembly.fasta |
| velvetg on data x_ Stats.tabular | HinfKW20_denovo_assembly_stats.txt |

6.3.0 Save your history as a workflow



Now that you've saved your files it's time to convert your history (which is also saved in Galaxy) into a workflow. A workflow is the equivalent of a bioinformatic pipeline – a series of commands that will run on your data to produce a final output.

Saving workflows will allow you to perform the same analyses that we just performed on a completely different dataset. For example, if there was a second *H. influenzae* strain that you had sequenced, you could run the same analysis we just ran on the *H. influenzae* KW20 strain with the same parameters in a single step. This is important for both efficiency and reproducibility.

6.3.1 Go to the [History pane](#) and select the [History options](#) button.

6.3.2 Under the [CURRENT HISTORY](#) section select

Notice that the Upload File tools include a description "This tool cannot be used in workflows". These will inevitably be removed from the workflow when we create it.

6.3.3 Simplify the [Workflow name](#) to **AssemblyandAnnotation**


A screenshot of the Galaxy web interface. The top navigation bar shows "Galaxy / /dguttman/galaxy547" and tabs for "Analyze Data", "Workflow", "Visualize", "Shared Data", "Help", "User", and a grid icon. On the left, a "Tools" sidebar lists categories like "GENOMICS ANALYSIS" and "METAGENOMICS". The main panel shows a "Workflow name" input field with the text "AssemblyandAnnotation". Below it are buttons for "Create Workflow", "Check all", and "Uncheck all". A "Tool" section lists "Upload File" (grayed out with the note "This tool cannot be used in workflows") and "FastQC" (checked with the note "Include 'FastQC' in workflow"). On the right, a "History Items created" table lists four items: "1 HinfKW20_rna.fna", "2 HinfKW20_datasubset.fastq", "3 FastQC on data 2: Webpage", and "4 FastQC on data 2: RawData". Each item has a checkbox and a "Treat as input dataset" field.

6.3.4 Select the [Create Workflow](#) button.


6.4.0 Editing your workflow

Now that a workflow has been created you can edit and rearrange portions as you see fit. You can even alter some of the parameters, add/remove steps, specify the files you want to include as output and resave the workflow or save it as a completely new one.

6.4.1 On the Galaxy menu select Workflows

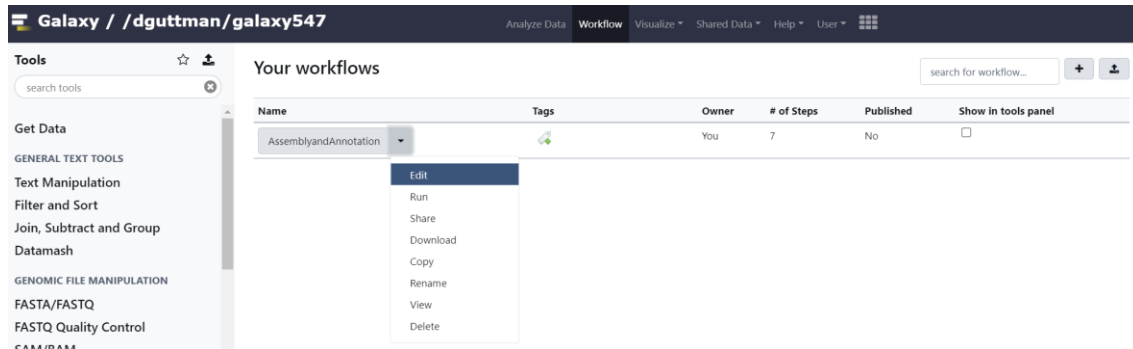

6.4.2 On the Workflows page select the pulldown menu for AssemblyandAnnotation and select the [redacted]. This will load the workflow editor as follows.

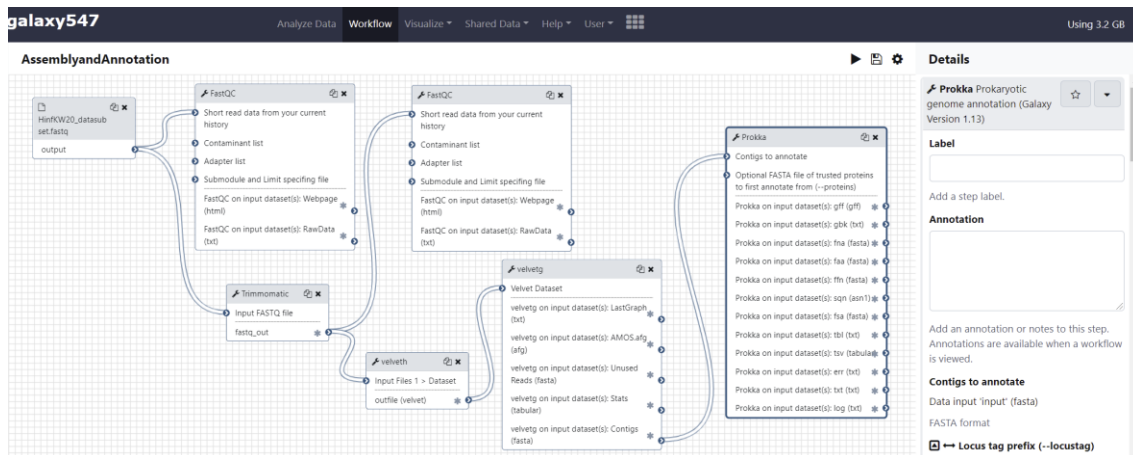

6.4.3 Use the visual editor to re-arrange the steps into a clear set of analyses. This will make it easier to see how you might want to modify the input/output files. Notice how the tools are connecting the output of one with the input of another.

6.4.4 In the right-hand Details pane you can see the specific parameters set for each tool when you click on the individual tool information boxes.

Notice that in this pane you can also label and annotate the individual steps. This can be a useful way to comment on what you are doing and why.

6.4.5 Click on the [redacted]

In the Details pane you can now see additional tags and annotations you can add to the overall workflow. These are helpful attributes to add if you are sharing your workflow with others as it will make it searchable and descriptive. You can also see previous saved versions of this workflow in the Version dropdown menu.

6.4.6 Use the Save icon (floppy disk) to save your workflow.

6.5.0 Run your saved workflow

Now that we've saved our updated workflow, how do we apply this to another dataset? Recall that our saved workflow appears to already include a dataset.

- 6.5.1 On the Galaxy menu select Workflows
- 6.5.2 On the Workflows page select the pulldown menu for AssemblyandAnnotation and select the Run option. This will load the workflow in it's run order as follows.
- 6.5.3 Under the History Options section, decide whether you'd like the results of this run to go to your current history or a new one. If you choose **Yes**, then be sure to give your history a name.
- 6.5.4 Step 1 will ask you to enter an initial file for analysis. The dropdown menu will include any available fastq files that are present in your current History pane. It is here that you will choose a set of fastq reads that are available in your history.

The screenshot shows the Galaxy web interface. The top navigation bar includes 'Galaxy / /duttman/galaxy547' and tabs for 'Analyze Data', 'Workflow', 'Visualize', 'Shared Data', 'Help', and 'User'. On the left, a 'Tools' sidebar lists categories like 'Get Data', 'GENERAL TEXT TOOLS', 'GENOMIC FILE MANIPULATION', and 'FASTA/FASTQ'. The main content area is titled 'Workflow: AssemblyandAnnotation' and features a 'Run workflow' button. Below this, the 'History Options' section allows users to choose whether to 'Send results to a new history' (Yes/No) and provides a field for 'History name' (currently 'AssemblyandAnnotation Session2'). A workflow diagram is visible, showing steps 1 through 5. Step 1 is 'HinfKW20_datasubset.fastq'. Step 2 is 'FastQC (G...'. Step 5 is 'Trimmomatic on HinfKW20_datasubset.fastq'. A dropdown menu for step 2 is open, showing '2: HinfKW20_datasubset.fastq' as the selected option. Below the workflow diagram, the 'Short read data' section shows 'Output dataset "output" from step 1'.

- 6.5.5 Double check the settings/parameters across your tool. You'll see that they can be edited individually but if you are saving this as a workflow, you'll want to keep them constant for consistency. If you're experimenting with your settings, you should note they will not be saved and you should experiment with these in the **Workflow Edit** page instead.

7.0.0 Class summary

That concludes our second lecture and introduction to working with the Galaxy platform. Next week we'll dig into some other tools on the platform to explore RNAseq. Altogether we've explored the following:

- Importing data files into Galaxy.
- Quality control analysis of raw and filtered reads with FastQC.
- Trimming and filter reads with Trimmomatic.
- Small genome assembly using Velvet.
- Prokaryotic genome annotation with Prokka.
- Downloading your datasets.
- Saving, editing, and running workflows based on your previous history.

4.1.0 Post-lecture assessment (10% of final grade)

Soon after this lecture, a homework assignment will be made available on Quercus in the assignment section. It will build on the ideas and/or data generated within this lecture. Each homework assignment will be worth 10% of your final mark. If you have assignment-related questions, please try the following steps in the order presented:

- Check the internet for a solution – read forums and learn to navigate for answers.
- Generate a discussion on Quercus outlining what you've tried so far and see if other students can contribute to a solution.
- Contact course teaching assistants or the instructor.

4.2.0 Suggested class preparation for Module 2

Next week we will begin exploring the Galaxy interface and learning how to perform quality control and genome assembly on data from sources like the SRA. To prepare for this, we suggest the following Coursera Modules:

- Genomic Data Science with Galaxy, James Taylor PhD:
 - Module 4: RNAseq (20 mins)

We also suggest the following reading which will be available from the Module 2 section of Quercus:

- *Altmann et al., 2012, Human Genetics, 131:1541-1554.*

